

BIOSAFETY MOBILE UNIT



The laboratory is a positive-pressure cleanroom system, with the operating room system protecting the operator and a design for storage examined until transferred to the laboratory.

**Overcome
restrictions
on field visits**

**Quickly
and
safely**

**Convenient for both
service providers
and recipients**

BIOSAFTEY MOBILE UNIT BMU-1STATION

FLUID SPRAY SYSTEM



ANTEROOM



GLOVE PORT



HEPA FILTER



SAMPLE STORAGE



INFECTIOUS
WASTE PASSAGES



Safety

- ✓ Positive pressure clean room control the pressure inside the chamber is always positive
- ✓ The air will flow from the top of booth and go through HEPA (H 13) filter to provide clean air environment inside



Confident

- ✓ Sample collection operations designed for operators to collect samples by inserting hands through rubber gloves
- ✓ Room disinfection system with ozone sterilizer



Secure

- ✓ A sample storage cabinet that maintains the temperature to the standards CDC
- ✓ CCTV system to check the surrounding situation while working
- ✓ Communicate with service recipients through the broadcasting system



Optional

- ✓ Internet WIFI systems within the car provide a variety of connection Application
- ✓ Install a Mobile power generator can work even without a power station

BIOSAFETY MOBILE UNIT BMU-2STATION

HEPA FILTER



GLOVE PORT




ANTEROOM




The pressure display screen in cleanroom


FLUID SPRAY SYSTEM

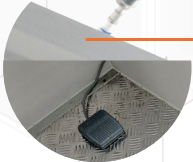

2 INFECTIOUS
WASTE PASSAGES

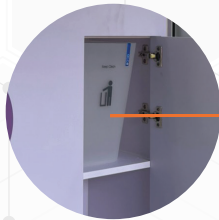

SAMPLE STORAGE

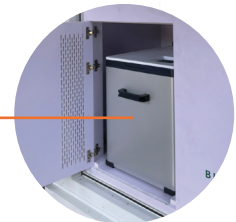

Safety

- ✓ Positive pressure clean room control the pressure inside the chamber is always positive
- ✓ The air will flow from the top of booth and go through HEPA (H 13) filter to provide clean air environment inside

Confident

- ✓ Sample collection operations designed for operators to collect samples by inserting hands through rubber gloves
- ✓ Room disinfection system with ozone sterilizer

Secure

- ✓ A sample storage cabinet that maintains the temperature to the standards CDC
- ✓ CCTV system to check the surrounding situation while working
- ✓ Communicate with service recipients through the broadcasting system

Optional

- ✓ Internet WIFI systems within the car provide a variety of connection Application
- ✓ Install a Mobile power generator can work even without a power station